

Eurosafe 60

Cleaning, disinfection
and drying
in complete safety

Eurosafe 60 pre-washes, washes, thermally disinfects and dries in a single cycle, so there is no need to clean by hand. This guarantees more effective and efficient results in less time and with fewer risks for the operator. The large, 60-litre tank can be used to thermally disinfect many items, including handpieces, loose instruments, cassettes and containers. The wide range of accessories optimizes the use of the space available. The user-friendly interface can be utilized to program up to 40 cycles, customizing them to match specific needs and materials by adjusting settings such as washing time and temperature. Eurosafe 60 guarantees complete, secure traceability and it becomes even more effective if connected to the other Euronda Pro System devices. It validates and records the data of all disinfection cycles, which can be saved on a USB key and transferred to a computer. Eurosafe 60 is available in two versions: with one or two peristaltic pumps.

**Classification:**

Medical Device Class IIb

**Dimensions (wxhxd):**

560 x 825 x 577 mm

**Weight:**

62 Kg

**Total power:**

2750 W 230 V

Pump power:

150 W

**Operating temperature:**

+5°C / +40°C



Support with wand holders



Injection washing system



Chemical liquid dispensing system



All in one

Eurosafe 60 pre-washes, washes, thermally disinfects and dries in a single cycle, for more effective results in less time.



Inspectionability

It is easy to check the contents during a cycle thanks to the tempered glass door.



Capacity

The large, 60-litre tank can thermally disinfect many items, including wands, loose instruments and containers.



Power

The injection washing system blasts water jets from a number of directions for outstanding cleaning results.



Eco-friendly

It only consumes 1.7 kW and 12 litres of water during a "Short" cycle.



Traceability

The disinfection process is validated and recorded using USB key, Ethernet connection, external printer.



Effectiveness

The forced air drying system removes all traces of water from hollow instruments, for optimum drying.



User-friendly

Up to 40 cycles can be programmed and password-protected thanks to the user-friendly interface of the LCD display.



Starter Kit



- 1 handbook
- 1 USB key
- 1 water test
- 1 Kg salt
- Euroclean 120 (500 ml)
- Eurobright 360 (500 ml, included just in the two-pump version)



Standard load basket



Insert for 8 trays
(40 mm slots)



Insert for vertically
positioned instruments



Basket for burs
and small parts



Eurosafe 60

Accessories



Load basket with 4 wand holders, 5 injectors for hollow items and 2 connections for aspiration tubes



Load basket with 6 wand holders and 6 injectors for hollow items



Load basket with 8 wand holders



Load basket with 11 injectors for hollow items and 2 connections for aspiration tubes



Net basket DIN 1/2



Net basket DIN 1/1




Net basket DIN 1/1 for rigid instruments



Turbine adaptor



Insert for 2 instrument containers or trays (50 mm slots)

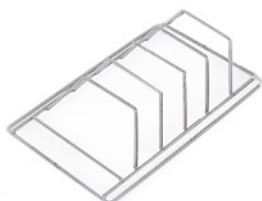

Insert for 4 instrument containers (35 mm slots)



External thermal paper printer



Adapters:

- male thread M3 X 0.5
- female thread M3 X 0.5
- female thread M3.5 X 0.6
- injector
- holder



Euroclean 120 - 3 L Enzymatic detergent for mechanic washing of medical devices



Eurobright 360 - 3 L Rinsing liquid, lubricant for medical equipment and devices in the washing machine



Datalogger
It transfers registered data from the thermo-desinfector to the PC, for a complete and guaranteed data traceability